



DESS

> DENTAL SMART SOLUTIONS



DESS

LIBRARY GUIDE

FOR HYBRID ABUTMENTS, BRIDGES, BARS AND
PROSTHESIS

Use of DESS Abutments Library in CAD dental systems

DESS Scan Abutments allow scans on implant level or abutment level. This depends on the type of connection and custom impression or intraoral scan.

The choice of the appropriate library connection depends also on the design of the prosthesis to be milled. DESS Library allows restorations on DESS Ti base, DESS AURUMBase® or restorations to be screwed direct on Straumann® synOcta® abutment, Nobel Biocare Multi- unit® abutment or ASTRA UniAbutment®.

Please note that in case of direct milling the connection is non-engaging only.

Depending on the CAD system, the library structure and connection names differ slightly from each other.

3Shape Dental System – implant level scan



For 3Shape Dental System on this level the following categories are available:

- DESS Ti base – for restorations connected with the implant via the Titanium Ti base. The scanning needs to be performed with the Scan Body for the specific implant connection. **Note:** This category includes parts **scanned from abutment level** - Multi-Unit abutment, UniAbutment® and synOcta® abutment - with its specific Scan Abutment (see Table 1.1).
- DESS MULTIUNIT & UNIABUTMENT PLUS Ti base – for restorations connected with the transepithelial abutment via the Ti base. The scanning needs to be performed with the Scan Body for the specific **implant connection**. The desired Multi-Unit abutment or UniAbutment® will be mounted virtually in the CAD Dental System. Please note that in this case the Ti base closing line (between the prosthesis and Ti base chamfer) depends on the proper kit selection (abutment height in mm),
- DESS AURUMBase® – for restorations connected with the implant via the AURUMBase®. The scanning needs to be performed with the Scan Body for the specific implant connection. **Note:** This category includes parts **scanned from abutment level** - Multi-Unit abutment - with its specific Scan Abutment (see Table 1.3).
- DESS Pre-milled Ti – for restorations connected with the implant via the pre-milled abutment made in Titanium. The scanning needs to be performed with the Scan Body for the specific implant connection. Note: The maximum width of the piece is restricted to 10 mm.
- DESS Pre-milled CoCr – for restorations connected with the implant via the pre-milled abutment made in CoCr. The scanning needs to be performed with the Scan Body for the specific implant connection. Note: The maximum width of the piece is restricted to 14 mm.

Take care to select the correct DESS Scan Abutment and the corresponding library file, as otherwise the milled prosthesis can differ from desired outcome and errors can occur, for example the level of Ti base settlement.

The tables on the following pages gives an overview of all DESS 3Shape Dental System library files.

Table 1.1 DESS TI-Base

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Nobel Active®	3.0	14.068		DESS TI-BASE	DESS NOB ACT TI-BASE	NOB_ACT_3.0 NE 15.068	52.068	
						53.100		
						NOB_ACT_3.0 NE 16.068	52.068	
						53.200		
	NP/3.5	14.041		DESS TI-BASE	DESS NOB ACT TI-BASE	NOB_ACT_NP NE 15.041 15.341	50.041	52.041
						53.100		
						NOB_ACT_NP E 16.041 16.341	50.041	52.041
						53.200		
	RP/4.3	14.042		DESS TI-BASE	DESS NOB ACT TI-BASE	NOB_ACT_RP NE 15.042 15.342	50.042	52.042
						53.100		
						NOB_ACT_RP E 16.042 16.342	50.042	52.042
						53.200		
WP/5.5	14.069		DESS TI-BASE	DESS NOB ACT TI-BASE	NOB_ACT_WP NE 15.069	52.069		
					53.100			
					NOB_ACT_WP NE 16.069	52.069		
					53.200			
Nobel Replace Select™	NP/3.5	14.004		DESS TI-BASE	DESS NOB RPL TI-BASE	NOB_RPL_NP NE 15.004 15.204 15.304	50.004	52.004
						53.100		
						NOB_RPL_NP E 16.004 16.204 16.304	50.004	52.004
						53.200		
	RP/4.3	14.005		DESS TI-BASE	DESS NOB RPL TI-BASE	NOB_RPL_RP NE 15.005 15.205 15.305	50.005	52.005
						53.100		
						NOB_RPL_RP E 16.005 16.205 16.305	50.005	52.005
						53.200		
	WP/5.0	14.006		DESS TI-BASE	DESS NOB RPL TI-BASE	NOB_RPL_WP NE 15.006 15.206 15.306	50.006	52.006
						53.100		
						NOB_RPL_WP E 16.006 16.206 16.306	50.006	52.006
						53.200		
6.0	14.097		DESS TI-BASE	DESS NOB RPL TI-BASE	NOB_RPL_NP 6.0 15.097	52.097		
					53.100			
					NOB_RPL_NP 6.0 16.097	52.097		
					53.200			
Nobel Brånemark®	NP/3.5	14.001		DESS TI-BASE	DESS NOB BRK TI-BASE	NOB_BRK_NP NE 15.001 15.201 15.301	50.001	52.001
						53.100		
						NOB_BRK_NP E 16.001 16.201 16.301	50.001	52.001
						53.200		
	RP/4.1	14.002		DESS TI-BASE	DESS NOB BRK TI-BASE	NOB_BRK_RP NE 15.002 15.202 15.302	50.002	52.002
						53.100		
						NOB_BRK_RP E 16.002 16.202 16.302	50.002	52.002
						53.200		
	WP/5.0	14.003		DESS TI-BASE	DESS NOB BRK TI-BASE	NOB_BRK_WP NE 15.003 15.203 15.303	50.003	
						52.003		
						NOB_BRK_WP E 16.003 16.203 16.303	50.003	
						52.003		

non-engaging

non-engaging

 Scan Abutment for desktop scanning
50.0XX

 Scan Abutment for intra-oral scanning
52.0XX

 Scan Abutment on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments			
3i Osseotite®	NP/3.4	14.011		DESS TI-BASE	DESS 3i OSS TI-BASE	3i_OSS_NP NE 15.011 15.211 15.311	50.011 52.011 53.100			
						3i_OSS_NP E 16.011 16.211 16.311	50.011 52.011 53.200			
						3i_OSS_RP NE 15.012 15.212 15.312	50.012 52.012 53.100			
						3i_OSS_NP E 16.011 16.212 16.312	50.012 52.012 53.200			
	RP/4.1	14.012		DESS TI-BASE	DESS 3i OSS TI-BASE	3i_OSS_WP NE 15.013 15.213 15.313	50.013 52.013 53.100			
						3i_OSS_WP E 16.013 16.213 16.313	50.013 52.013 53.200			
	3i Certain®	NP/3.4	14.014		DESS TI-BASE	DESS 3i CER TI-BASE	3i_CER_NP NE 15.014 15.214 15.314	50.014 52.014 53.100		
							3i_CER_NP E 16.014 16.214 16.314	50.014 52.014 53.200		
							3i_CER_RP NE 15.015 15.215 15.315	50.015 52.015 53.100		
							3i_CER_RP E 16.015 16.215 16.315	50.015 52.015 53.200		
RP/4.1		14.015		DESS TI-BASE	DESS 3i CER TI-BASE	3i_CER_WP NE 15.016 15.216 15.316	50.016 52.016 53.100			
						3i_CER_WP E 16.016 16.216 16.316	50.016 52.016 53.200			
Straumann® Soft Tissue Level		NNC/3.5	14.070		DESS TI-BASE	DESS STR STLV TI-BASE	NOB_ACT_WP NE 15.069	52.069 53.100		
							NOB_ACT_WP NE 16.069	52.069 53.200		
							STR STLV RN NE 15.046	50.009 52.009		
							STR STLV RN E 16.009	50.009 52.009		
	RN on synOCTA®	14.010		DESS TI-BASE	DESS STR STLV TI-BASE	STR SYN RN NE 15.010	50.010 52.010			
						STR SYN RN E 16.010	50.010 52.010			
	WN/6.5	14.045		DESS TI-BASE	DESS STR STLV TI-BASE	STR STLV WP NE 15.045	50.045 52.045			
						STR STLV WP E 16.045	50.045 52.045			
WN on synOCTA®	14.046		DESS TI-BASE	DESS STR STLV TI-BASE	STR SYN WN NE 15.048	50.015 52.015				
					STR SYN WN E 16.048	50.015 52.015				

non-engaging

non-engaging

 Scan Abutment
for desktop scanning
50.0XX

 Scan Abutment
for intra-oral scanning
52.0XX

 Scan Abutment
on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Straumann® Bone Level	NC/3.3	14.043		DESS TI-BASE	DESS STR BLV TI-BASE	STR_BLV_NC NE 15.043 15.343	50.043	52.043
						TR_BLV_NC E 16.043 16.343	50.043	52.043
	RC/4.1	14.044		DESS TI-BASE	DESS STR BLV TI-BASE	STR_BLV_RC NE 15.044 15.344	50.044	52.044
						TR_BLV_RC E 16.044 16.344	50.044	52.044
	SRA	14.098		DESS TI-BASE	DESS STR BLV TI-BASE	STR_BLV_SRA 15.098	52.098	
						STR_BLV_SRA 16.098	52.098	
Zimmer Screw-Vent®	NP/3.5	14.017		DESS TI-BASE	DESS ZIM SCV TI-BASE	ZIM_SCV_NP NE 15.017 15.217 15.317	50.017	52.017
						ZIM_SCV_NP E 16.017 16.217 16.317	50.017	52.017
	RP/4.5	14.018		DESS TI-BASE	DESS ZIM SCV TI-BASE	ZIM_SCV_RP NE 15.018 15.218 15.318	50.018	52.018
						ZIM_SCV_RP E 16.018 16.218 16.318	50.018	52.018
	WP/5.7	14.019		DESS TI-BASE	DESS ZIM SCV TI-BASE	ZIM_SCV_WP NE 15.019 15.219 15.319	50.019	52.019
						ZIM_SCV_WP E 16.019 16.219 16.319	50.019	52.019
Astra tech implant system™ EV	EV/3.6	14.059		DESS TI-BASE	DESS AST EVO TI-BASE	AST_EVO 3.6_NE 15.059 15.359	50.059	52.059
						AST_EVO 3.6_E 16.059 16.359	50.059	52.059
	EV/4.2	14.060		DESS TI-BASE	DESS AST EVO TI-BASE	AST_EVO 4.2_NE 15.060 15.360	50.060	52.060
						AST_EVO 4.2_E 16.060 16.360	50.060	52.060
	EV/4.8	14.061		DESS TI-BASE	DESS AST EVO TI-BASE	AST_EVO 4.8_NE 15.061 15.361	50.061	52.061
						AST_EVO 4.8_E 16.061 16.361	50.061	52.061
	Uniabutment®	14.058		DESS TI-BASE	DESS AST EVO TI-BASE	AST_EVO uniabutment_NE 15.058	50.058	52.058
						AST_EVO uniabutment_NE 16.058	50.058	52.058

 non-engaging

 non-engaging

 Scan Abutment
for desktop scanning
50.0XX

 Scan Abutment
for intra-oral scanning
52.0XX

 Scan Abutment
on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Astra tech osseospeed™	NP/3.0	14.023		DESS TI-BASE	DESS AST OSS TI-BASE	AST_OSS_NP NE 15.023	52.023	
						AST_OSS_NP E 16.023	50.059	52.059
	RP/3.5-4.0	14.024		DESS TI-BASE	DESS AST OSS TI-BASE	AST_OSS_RP NE 15.024 15.324	50.024	52.024
						AST_OSS_RP E 16.024 16.324	50.024	52.024
	WP/4.5-5.0	14.025		DESS TI-BASE	DESS AST OSS TI-BASE	AST_OSS_WP NE 15.025 15.325	50.025	52.025
						AST_OSS_WP E 16.025 16.325	50.025	52.025
	Uniabutment® 20°	14.036		DESS TI-BASE	DESS AST OSS TI-BASE	AST_OSS_Uniabutment®_20° NE 15.036	50.036	52.036
Friadent XIVE®	Uniabutment® 45°	14.037		DESS TI-BASE	DESS AST OSS TI-BASE	AST_OSS_Uniabutment®_45° NE 15.037	50.037	52.037
	NP/3.4	14.038		DESS TI-BASE	DESS DEN XIV TI-BASE	DEN_XIV_NP NE 15.038	50.038	52.038
						DEN_XIV_NP E 16.038	50.038	52.038
	RP/3.8	14.039		DESS TI-BASE	DESS DEN XIV TI-BASE	DEN_XIV_RP NE 15.039	50.039	52.039
						DEN_XIV_RP E 16.039	50.039	52.039
	WP/4.5	14.040		DESS TI-BASE	DESS DEN XIV TI-BASE	DEN_XIV_WP NE 15.040	50.040	52.040
						DEN_XIV_WP E 16.040	50.040	52.040
Ankylos® C/X	All platforms, height GH 1,5	14.048 14.049		DESS TI-BASE	DESS DEN ANK TI-BASE	DEN ANK NE GH1 15.050	50.049	52.049
						DEN ANK E GH1 16.050	50.049	52.049
	All platforms, height GH 2,5	14.048 14.049		DESS TI-BASE	DESS DEN ANK TI-BASE	DEN ANK NE GH2 15.250	50.049	52.049
						DEN ANK E GH2 16.250	50.049	52.049

 non-engaging

 non-engaging

 Scan Abutment for desktop scanning
50.0XX

 Scan Abutment for intra-oral scanning
52.0XX

 Scan Abutment on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Osstem® TS	Mini Yellow	14.072		DESS TI-BASE	DESS OSS TS TI-BASE	OSS_TS MINI_NE 15.072	50.072	52.072
						15.372	53.100	
	STD Regular	14.073		DESS TI-BASE	DESS OSS TS TI-BASE	OSS_TS MINI_E 16.072	50.072	52.072
						16.372	53.200	
Camlog®	3.3	14.064		DESS TI-BASE	DESS CAMLOG TI-BASE	CAM_3.3 NE 15.064	52.064	
							53.100	
						CAM_3.3 E 16.064	52.064	
							53.200	
	3.8	14.065		DESS TI-BASE	DESS CAMLOG TI-BASE	15.065	50.065	52.065
						CAM_3.8 NE 15.265	53.100	
						15.365		
						CAM_3.8 E 16.065	50.065	52.065
						16.265	53.200	
						16.365		
	4.3	14.066		DESS TI-BASE	DESS CAMLOG TI-BASE	15.066	50.066	52.066
						CAM_4.3 NE 15.266	53.100	
Conelog®	NP/3.3	14.074		DESS TI-BASE	DESS CONELOG TI-BASE	15.366		
						CAM_4.3 E 16.066	50.066	52.066
						16.066	53.200	
						16.266		
						16.366		
	5.0	14.067		DESS TI-BASE	DESS CAMLOG TI-BASE	CAM_5.0 NE 15.067	52.067	
							53.100	
						CAM_5.03 E 16.067	52.067	
							53.200	
	RP/3.8	14.075		DESS TI-BASE	DESS CONELOG TI-BASE	CONELOG NP NE 15.074	52.074	
							53.100	
Conelog®	WP/5.0	14.076		DESS TI-BASE	DESS CONELOG TI-BASE	CONELOG NP E 16.074	52.074	
							53.200	
						CONELOG RP NE 15.075	52.075	
							53.100	
						CONELOG RP E 16.075	52.075	
							53.200	
						CONELOG WP NE 15.076	52.076	
							53.100	
						CONELOG WP E 16.076	52.076	
							53.200	

 non-engaging

 non-engaging

 Scan Abutment for desktop scanning
50.0XX

 Scan Abutment for intra-oral scanning
52.0XX

 Scan Abutment on TI base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Biohorizns® Internal	NP/3.5	14.087	 	DESS TI-BASE	DESS BH INT TI-BASE	BH_INTERNAL_NP NE 15.087	52.087	
						BH_INTERNAL_NP E 16.087	53.100	
	RP/4.5	14.088	 	DESS TI-BASE	DESS BH INT TI-BASE	BH_INTERNAL_RP NE 15.088	52.088	
						BH_INTERNAL_RP E 16.088	53.200	
	WP/5.7	14.089	 	DESS TI-BASE	DESS BH INT TI-BASE	BH_INTERNAL_WP NE 15.089	52.089	
						BH_INTERNAL_WP E 16.089	53.200	
Biohorizns® External	NP/3.7	14.020	 	DESS TI-BASE	DESS BH EXT TI-BASE	BH_INTERNAL_NP NE 15.051	50.020	52.020
						BH_INTERNAL_NP E 16.051	53.100	
	RP/4.2	14.021	 	DESS TI-BASE	DESS BH EXT TI-BASE	BH_EXTERNAL_RP NE 15.052	50.021	52.021
						BH_EXTERNAL_RP E 16.052	53.200	
	WP/5.2	14.022	 	DESS TI-BASE	DESS BH EXT TI-BASE	BH_EXTERNAL_WP NE 15.053	50.022	52.022
						BH_EXTERNAL_WP E 16.053	53.200	
MIS® Seven	NP	14.056	 	DESS TI-BASE	DESS MIS SEV TI-BASE	MIS_SEV_NP_NE 15.056	50.056	52.056
						MIS_SEV_NP_E 16.056	53.100	
	SP	14.095	 	DESS TI-BASE	DESS MIS SEV TI-BASE	MIS_SEV_SP_NE 15.095 15.295 15.395	50.017	52.017
						MIS_SEV_SP_E 16.095 16.295 16.395	53.200	
	WP	14.096	 	DESS TI-BASE	DESS MIS SEV TI-BASE	MIS_SEV_WP_NE 15.096 15.296 15.396	50.018	52.018
						MIS_SEV_WP_E 16.096 16.296 16.396	53.200	

 non-engaging

 non-engaging

 Scan Abutment
for desktop scanning
50.0XX

 Scan Abutment
for intra-oral scanning
52.0XX

 Scan Abutment
on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments
MIS® C1 internal	NP	14.081		DESS TI-BASE	DESS MIS C1 TI-BASE	MIS C1 NP NE 15.081	52.081 53.100
						MIS C1 NP E 16.081	52.081 53.200
	SP	14.082		DESS TI-BASE	DESS MIS C1 TI-BASE	MIS C1 SP NE 15.082	52.082 53.100
						MIS C1 SP E 16.082	52.082 53.200
	WP	14.083		DESS TI-BASE	DESS MIS C1 TI-BASE	MIS C1 WP NE 15.083	52.083 53.100
						MIS C1 WP E 16.083	52.083 53.200
	Multiunit 4.8	14.085		DESS TI-BASE	DESS MIS C1 TI-BASE	MIS C1 MULTIUNIT WE 15.083	52.085
MIS® V3	NP	14.084		DESS TI-BASE	DESS MIS V3 TI-BASE	MIS C1 NP NE 15.084	52.084 53.100
						MIS C1 NP E 16.084	52.084 53.200
	SP	14.082		DESS TI-BASE	DESS MIS V3 TI-BASE	MIS C1 SP NE 15.082	52.083 53.100
						MIS C1 SP E 16.082	52.082 53.200
	Multiunit 4.8	14.085		DESS TI-BASE	DESS MIS V3 TI-BASE	MIS C1 MULTIUNIT 4.8 15.085	52.085
Megagen AntRidge®		14.057		DESS TI-BASE	DESS MGN ANY TI-BASE	MGN ANY NE 15.057	50.057 53.100
						MGN ANY E 16.057	50.057 53.200
Anthogyr AXIOM®		14.077		DESS TI-BASE	DESS ANTHOGRYR TI-BASE	ANTHOGYR Ti base NE 15.077	50.077 53.100
						ANTHOGYR Ti base E 16.077	50.077 53.200
Medentis ICX®		14.080		DESS TI-BASE	DESS MED ICX TI-BASE	MED ICX NE 15.080	50.080 53.100
						MED ICX E 15.080	50.080 53.200

 non-engaging

 non-engaging

 Scan Abutment
for desktop scanning
50.0XX

 Scan Abutment
for intra-oral scanning
52.0XX

 Scan Abutment
on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments
GlobalID IN-KONE®		14.086		DESS TI-BASE	DESS GLOBAL ID TI-BASE	IN_KONE Ti base NE 15.086	52.086
							53.100
						IN_KONE Ti base E 16.086	52.086
							53.200
Dentium Superline™ & Implantium		14.090		DESS TI-BASE	DESS DENTIUM TI-BASE	DENTIUM Ti base NE 15.090	52.090
							53.100
						DENTIUM Ti base E 15.090	52.090
							53.200
Bego Semados®	3.25/3.75	14.091		DESS TI-BASE	DESS BEGO TI-BASE	BEGO Ti base 3.25 NE 15.091	52.091
							53.100
						BEGO Ti base 3.25 E 16.091	52.091
							53.200
	4.1	14.092		DESS TI-BASE	DESS BEGO TI-BASE	BEGO Ti base 4.1 NE 15.091	52.092
							53.100
						BEGO Ti base 4.1 E 16.091	52.092
							53.200
	4.5	14.093		DESS TI-BASE	DESS BEGO TI-BASE	BEGO Ti base 4.5 NE 15.091	52.093
							53.100
						BEGO Ti base 4.5 E 16.091	52.093
							53.200
	5.5	14.094		DESS TI-BASE	DESS BEGO TI-BASE	BEGO Ti base 5.5 NE 15.094	52.094
							53.100
						BEGO Ti base 5.5 E 16.094	52.094
							53.200
Biotech Kontakt®		14.078		DESS TI-BASE	DESS BIOTECH TI-BASE	BIOTECH Ti base NE 15.078	52.078
							53.100
						BIOTECH Ti base E 16.078	52.078
							53.200
BTI®		14.079		DESS TI-BASE	DESS BTI TI-BASE	BTI INTERNAL Base NE 15.079	52.079
							53.100
						BTI INTERNAL Base E 15.079	52.079
							53.200
Multi-Unit	RP	14.007		DESS TI-BASE	DESS Multi-Unit TI-BASE	Multi_Unit_RP NE 15.007	50.007
							52.007
	RP	14.007		DESS TI-BASE	DESS Multi-Unit TI-BASE	Multi_Unit_RP NE 15.007/8	50.007
							52.007
	WP	14.047		DESS TI-BASE	DESS Multi-Unit TI-BASE	Multi_Unit_RP E 15.054	50.047
							52.047

 non-engaging

 non-engaging

 Scan Abutment
for desktop scanning
50.0XX

 Scan Abutment
for intra-oral scanning
52.0XX

 Scan Abutment
on Ti base
53.100 53.200

**Table 1.2 DESS
scan on Uniabutment
& Multi-Unit &
Synocta**

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments
Uniabutment® on Astra Tech™ Osseospeed™	Uniabutment® 20°	14.036	●	DESS UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	DESS AST UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	20	50.036 52.036
	Uniabutment® 45°	14.037	●	DESS UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	DESS AST UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	45	50.037 52.037
Astra Tech™ Implant System™ EV	Uniabutment®	14.058	●	DESS UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	DESS AST UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	EV Uniabutment	50.058 52.058

Multi-Unit	RP	14.007	●	DESS MULTIUNIT DIRECT (Multi-Unit Scan Abutment)	DESS AST UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	Multi-Unit	50.007 52.007
	WP	14.047	●	DESS MULTIUNIT DIRECT (Multi-Unit Scan Abutment)	DESS AST UNIABUTMENT DIRECT (Uniabutment Scan Abutment)	Multi-Unit WP	50.047 52.047

Straumann® Soft Tissue level witch synOcta® Abutment	RN on synOcta®	14.010	●	DESS STR SYN DIRECT (Bars and bridges)	DESS STR SYN DIRECT (Bars and bridges)	STR SYN RN 51.010	50.010 52.010
	WN on synOcta®	14.046	●	DESS STR SYN DIRECT (Bars and bridges)	DESS STR SYN DIRECT (Bars and bridges)	STR SYN WN 51.048	50.046 52.037

● non-engaging

● non-engaging

Scan Abutment
for desktop scanning
50.OXX

Scan Abutment
for intra-oral scanning
52.OXX

Table 1.3 DESS ANGLEBase®

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Nobel Active®	NP	14.041		DESS AURUMBase	DESS NOB ACT AURUMBase	NOB_ACT_NP NE 36.041	50.041	52.041
							53.100	
	RP	14.042		DESS AURUMBase	DESS NOB ACT AURUMBase	NOB_ACT_NP E 37.041	50.041	52.041
							53.200	
Nobel Replace Select™	NP	14.004		DESS AURUMBase	DESS NOB RPL AURUMBase	NOB_ACT_RP NE 36.042	50.042	52.042
							53.100	
	RP	14.005		DESS AURUMBase	DESS NOB RPL AURUMBase	NOB_ACT_RP E 37.042	50.042	52.042
							53.200	
	WP	14.006		DESS AURUMBase	DESS NOB RPL AURUMBase	NOB_RPL_NP NE 36.004	50.004	52.004
							53.100	
Nobel Brånemark®	NP	14.001		DESS AURUMBase	DESS NOB BRK AURUMBase	NOB_RPL_NP E 37.004	50.004	52.004
							53.200	
	RP	14.002		DESS AURUMBase	DESS NOB BRK AURUMBase	NOB_RPL_RP NE 36.005	50.005	52.005
							53.100	
	WP	14.003		DESS AURUMBase	DESS NOB BRK AURUMBase	NOB_RPL_RP E 37.005	50.005	52.005
							53.200	
	NP	14.001		DESS AURUMBase	DESS NOB BRK AURUMBase	NOB_RPL_WP NE 36.006	50.006	52.006
							53.100	
3i Osseotite®	NP/3.4	14.011		DESS AURUMBase	DESS NOB BRK AURUMBase	NOB_RPL_WP E 37.006	50.006	52.006
							53.200	
	RP/4.1	14.012		DESS AURUMBase	DESS 3i OSS AURUMBase	NOB_BRK_NP NE 36.001	50.001	52.001
							53.100	
	WP/5.0	14.013		DESS AURUMBase	DESS 3i OSS AURUMBase	NOB_BRK_NP E 37.001	50.001	52.001
							53.200	
	NP/3.4	14.011		DESS AURUMBase	DESS 3i OSS AURUMBase	NOB_BRK_RP NE 36.002	50.002	52.002
							53.100	
	RP/4.1	14.012		DESS AURUMBase	DESS 3i OSS AURUMBase	NOB_BRK_RP E 37.002	50.002	52.002
							53.200	
3i Osseotite®	NP/3.4	14.011		DESS AURUMBase	DESS 3i OSS AURUMBase	NOB_BRK_WP NE 36.003	50.003	52.003
							53.100	
	RP/4.1	14.012		DESS AURUMBase	DESS 3i OSS AURUMBase	NOB_BRK_WP E 37.003	50.003	52.003
							53.200	
	WP/5.0	14.013		DESS AURUMBase	DESS 3i OSS AURUMBase	3I_OSS_NP NE 36.011	50.011	52.011
							53.100	
3i Osseotite®	NP/3.4	14.011		DESS AURUMBase	DESS 3i OSS AURUMBase	3I_OSS_NP E 37.011	50.011	52.011
							53.200	
	RP/4.1	14.012		DESS AURUMBase	DESS 3i OSS AURUMBase	3I_OSS_RP NE 36.012	50.012	52.012
							53.100	
	WP/5.0	14.013		DESS AURUMBase	DESS 3i OSS AURUMBase	3I_OSS_RP E 37.012	50.012	52.012
							53.200	
3i Osseotite®	NP/3.4	14.011		DESS AURUMBase	DESS 3i OSS AURUMBase	3I_OSS_WP NE 36.013	50.013	52.013
							53.100	
3i Osseotite®	RP/4.1	14.012		DESS AURUMBase	DESS 3i OSS AURUMBase	3I_OSS_WP E 37.013	50.013	52.013
							53.200	

 non-engaging

 non-engaging

 Scan Abutment for desktop scanning
50.0XX

 Scan Abutment for intra-oral scanning
52.0XX

 Scan Abutment on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
3i Certain®	NP/3.4	14.014		DESS AURUMBase	DESS 3i CER AURUMBase	3I_CER_NP NE 36.014	50.014	52.014
						3I_CER_NP E 37.014	53.100	
	RP/4.1	14.015		DESS AURUMBase	DESS 3i CER AURUMBase	3I_CER_RP NE 36.015	50.014	52.014
						3I_CER_RP E 37.015	53.200	
	WP/5.0	14.016		DESS AURUMBase	DESS 3i CER AURUMBase	3I_CER_WP NE 36.016	50.015	52.015
						3I_CER_WP E 37.016	53.100	
Straumann® Soft Tissue Level	RN/4.8	14.009		DESS AURUMBase	DESS STR STL AURUMBase	STR_STLV_RN NE 36.046	50.016	52.016
						STR_STLV_RN E 37.046	53.200	
	WN/6.5	14.045		DESS AURUMBase	DESS STR STL AURUMBase	STR_STLV_WP NE 36.045	50.015	52.015
						STR_STLV_WP E 37.045	53.100	
Straumann® Bone Level	NC/3.3	14.043		DESS AURUMBase	DESS STR BLV AURUMBase	STR_BLV_NC NE 36.043	50.009	52.009
						TR_BLV_NC E 37.043	53.100	
	RC/4.1	14.044		DESS TI-BASE	DESS STR BLV AURUMBase	STR_BLV_RC NE 36.044	50.009	52.009
						TR_BLV_RC E 37.044	53.200	
Zimmer Screw-Vent®	NP/3.5	14.017		DESS AURUMBase	DESS ZIM SCV AURUMBase	ZIM_SCV_NP NE 36.017	50.043	52.043
						ZIM_SCV_NP E 37.017	53.100	
	RP/4.5	14.018		DESS AURUMBase	DESS ZIM SCV AURUMBase	ZIM_SCV_RP NE 36.018	50.043	52.043
						ZIM_SCV_RP E 37.018	53.200	
	WP/5.7	14.019		DESS AURUMBase	DESS ZIM SCV AURUMBase	ZIM_SCV_WP NE 36.019	50.044	52.044
						ZIM_SCV_WP E 37.019	53.100	

 non-engaging

 non-engaging

 Scan Abutment
for desktop scanning
50.0XX

 Scan Abutment
for intra-oral scanning
52.0XX

 Scan Abutment
on Ti base
53.100 53.200

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Astra tech implant system™ EV	EV/3.6	14.059		DESS AURUMBase	DESS AST EVO AURUMBase	AST_EVO 3.6_NE 36.059	50.059	52.059
						AST_EVO 3.6_E 37.059	50.059	52.059
	EV/4.2	14.060		DESS AURUMBase	DESS AST EVO AURUMBase	AST_EVO 4.2_NE 36.060	50.060	52.060
						AST_EVO 4.2_E 37.060	50.060	52.060
	EV/4.8	14.061		DESS AURUMBase	DESS AST EVO AURUMBase	AST_EVO 4.8_NE 36.061	50.061	52.061
						AST_EVO 4.8_E 37.061	50.061	52.061
Astra tech osseospeed™	RP	14.024		DESS AURUMBase	DESS AST OSS AURUMBase	AST_OSS_RP NE 36.024	50.024	52.024
						AST_OSS_RP E 37.024	50.024	52.024
	WP	14.025		DESS AURUMBase	DESS AST OSS AURUMBase	AST_OSS_WP NE 36.025	50.025	52.025
						AST_OSS_WP E 37.025	50.025	52.025
Friadent XIVE®	NP	14.038		DESS AURUMBase	DESS DEN XIV AURUMBase	DEN_XIV_NP NE 36.038	50.038	52.038
						DEN_XIV_NP E 37.038	50.038	52.038
	RP	14.039		DESS AURUMBase	DESS DEN XIV AURUMBase	DEN_XIV_RP NE 36.039	50.039	52.039
						DEN_XIV_RP E 37.039	50.039	52.039
	WP	14.040		DESS AURUMBase	DESS DEN XIV AURUMBase	DEN_XIV_WP NE 36.040	50.040	52.040
						DEN_XIV_WP E 37.040	50.040	52.040
 non-engaging		 non-engaging		 Scan Abutment for desktop scanning 50.0XX	 Scan Abutment for intra-oral scanning 52.0XX	 Scan Abutment on Ti base 53.100 53.200		

	Platform	Analogue	Type	Category	System	Kit	Scan Abutments	
Osstem® TS	MINI	14.072		DESS AURUMBase	DESS OSS TS AURUMBase	OSS_TS MINI_NE 36.072	50.072	52.072
							53.100	
	STD	14.073		DESS AURUMBase	DESS OSS TS AURUMBase	OSS_TS MINI_E 37.072	50.072	52.072
							53.200	
Camlog®	3.8	14.065		DESS AURUMBase	DESS CAMLOG AURUMBase	CAM_3.8 NE 36.065	50.065	52.065
							53.100	
	4.3	14.066		DESS AURUMBase	DESS CAMLOG AURUMBase	CAM_3.8 E 36.065	50.065	52.065
							53.200	
Mis® Seven	SP	14.095		DESS AURUMBase	DESS MIS SEV AURUMBase	CAM_4.3 NE 37.066	50.066	52.066
							53.100	
	WP	14.096		DESS AURUMBase	DESS MIS SEV AURUMBase	CAM_4.3 E 37.066	50.066	52.066
							53.200	
Megagen AnyRidge®	SP	14.095		DESS AURUMBase	DESS MIS SEV AURUMBase	MIS_SEV_SP NE 36.017	50.017	52.017
							53.100	
	WP	14.096		DESS AURUMBase	DESS MIS SEV AURUMBase	MIS_SEV_SP E 37.017	50.017	52.017
							53.200	
Medentis ICX®	RP	14.007		DESS AURUMBase	DESS Multi-Unit AURUMBase	STR_BLV_RC WP NE 36.018	50.018	52.018
							53.100	
	WP	14.047		DESS AURUMBase	DESS Multi-Unit AURUMBase	STR_BLV_RC WP E 37.018	50.018	52.018
							53.200	
Multi-Unit	RP	14.007		DESS AURUMBase	DESS MGN ANY AURUMBase	MGN_ANY_NE 36.057	50.017	52.017
							53.100	
	WP	14.047		DESS AURUMBase	DESS MGN ANY AURUMBase	MGN_ANY_E 37.057	50.017	52.017
							53.200	
Medentis ICX®	RP	14.007		DESS AURUMBase	DESS MED ICX AURUMBase	MED_ICX_NE 36.080	50.080	52.080
							53.100	
	WP	14.047		DESS AURUMBase	DESS MED ICX AURUMBase	MED_ICX_E 37.080	50.080	52.080
							53.200	

 non-engaging

 non-engaging




 Scan Abutment for desktop scanning
50.0XX

 Scan Abutment for intra-oral scanning
52.0XX

 Scan Abutment on Ti base
53.100 53.200

**Table 1.4 DESS
Pre-milled
Abutment Ti**

	Platform	Analogue	Category	System	Kit	Scan Abutments	
Nobel Active®	3.0	14.068	DESS PRE-MILLED ABUTMENT TI	DESS NOB ACT PRE-MILLED TI	NOB_ACT_3.0 E 61.068	52.068	
	NP	14.041	DESS PRE-MILLED ABUTMENT TI	DESS NOB ACT PRE-MILLED TI	NOB_ACT_NP E 61.041 _{ø10mm} 62.041 _{ø14mm}	50.041	52.041
	RP	14.042	DESS PRE-MILLED ABUTMENT TI	DESS NOB ACT PRE-MILLED TI	NOB_ACT_RP E 61.042 _{ø10mm} 62.042 _{ø14mm}	50.042	52.042
	WP	14.069	DESS PRE-MILLED ABUTMENT TI	DESS NOB ACT PRE-MILLED TI	NOB_ACT_WP E 61.069	52.069	
Nobel Replace Select™	NP	14.004	DESS PRE-MILLED ABUTMENT TI	DESS NOB RPL PRE-MILLED TI	NOB_RPL_NP E 61.004	50.004	52.004
	RP	14.005	DESS PRE-MILLED ABUTMENT TI	DESS NOB RPL PRE-MILLED TI	NOB_RPL_RP E 61.005	50.005	52.005
	WP	14.006	DESS PRE-MILLED ABUTMENT TI	DESS NOB RPL PRE-MILLED TI	NOB_RPL_WP E 61.006	50.006	52.006
	6.0	14.097	DESS PRE-MILLED ABUTMENT TI	DESS NOB RPL PRE-MILLED TI	NOB_RPL_NP 6.0 E 61.097	52.097	
Nobel Brånemark®	NP	14.001	DESS PRE-MILLED ABUTMENT TI	DESS NOB BRK PRE-MILLED TI	NOB_BRK_NP E 61.001	50.001	52.001
	RP	14.002	DESS PRE-MILLED ABUTMENT TI	DESS NOB BRK PRE-MILLED TI	NOB_BRK_RP E 61.002	50.002	52.002
	WP	14.003	DESS PRE-MILLED ABUTMENT TI	DESS NOB BRK PRE-MILLED TI	NOB_BRK_WP E 61.003	50.003 52.003	
3i Osseotite®	NP	14.011	DESS PRE-MILLED ABUTMENT TI	DESS 3i OSS PRE-MILLED TI	3i_OSS_NP E 61.011	50.011	52.011
	RP	14.012	DESS PRE-MILLED ABUTMENT TI	DESS 3i OSS PRE-MILLED TI	3i_OSS_RP E 61.012	50.012	52.012
	WP	14.013	DESS PRE-MILLED ABUTMENT TI	DESS 3i OSS PRE-MILLED TI	3i_OSS_WP E 61.013	50.013	52.013
3i Certain®	NP	14.014	DESS PRE-MILLED ABUTMENT TI	DESS 3i CER PRE-MILLED TI	3i_CER_NP E 61.014	50.014	52.014
	RP	14.015	DESS PRE-MILLED ABUTMENT TI	DESS 3i CER PRE-MILLED TI	3i_CER_RP E 61.015	50.015	52.015
	WP	14.016	DESS PRE-MILLED ABUTMENT TI	DESS 3i CER PRE-MILLED TI	3i_CER_WP E 61.016	50.016	52.016
Straumann® Soft Tissue Level	NNC	14.070	DESS PRE-MILLED ABUTMENT TI	DESS STR STL PRE-MILLED TI	NOB_ACT_WP E 61.070	52.069	
	RN	14.009	DESS PRE-MILLED ABUTMENT TI	DESS STR STL PRE-MILLED TI	STR_STLV_RN E 61.009 _{ø10mm} 62.009 _{ø14mm}	50.009 52.009	
	WN	14.045	DESS PRE-MILLED ABUTMENT TI	DESS STR STL PRE-MILLED TI	STR_STLV_WN E 61.045 _{ø10mm} 62.045 _{ø14mm}	50.045 52.045	

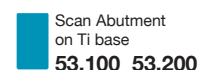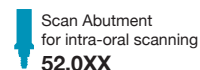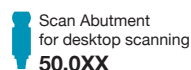

	Platform	Analogue	Category	System	Kit	Scan Abutments
Straumann® Bone Level	NC	14.043	DESS PRE-MILLED ABUTMENT TI	DESS STR BLV PRE-MILLED TI	STR_BLV_NC E 61.043 62.043	50.043 52.043
	RC	14.044	DESS PRE-MILLED ABUTMENT TI	DESS STR BLV PRE-MILLED TI	STR_BLV_RC E 61.044 62.044	50.044 52.044
Zimmer Screw-Vent®	NP	14.017	DESS PRE-MILLED ABUTMENT TI	DESS ZIM SCV PRE-MILLED TI	ZIM_SCV_NP E 61.017 62.017	50.017 52.017
	RP	14.018	DESS PRE-MILLED ABUTMENT TI	DESS ZIM SCV PRE-MILLED TI	ZIM_SCV_RP E 61.018	50.018 52.018
	WP	14.019	DESS PRE-MILLED ABUTMENT TI	DESS ZIM SCV PRE-MILLED TI	ZIM_SCV_WP E 61.019	50.019 52.019
Astra tech implant system™ EV	EV/3.6	14.059	DESS PRE-MILLED ABUTMENT TI	DESS AST EVO PRE-MILLED TI	AST_EVO 3.6_E 61.059	50.059 52.059
	EV/4.2	14.060	DESS PRE-MILLED ABUTMENT TI	DESS AST EVO PRE-MILLED TI	AST_EVO 4.2_E 61.060	50.060 52.060
	EV/4.8	14.061	DESS PRE-MILLED ABUTMENT TI	DESS AST EVO PRE-MILLED TI	AST_EVO 4.8_E 61.061	50.061 52.061
Astra tech implant Osseospeed™	RP	14.024	DESS PRE-MILLED ABUTMENT TI	DESS AST OSS PRE-MILLED TI	AST_OSS_RP E 61.024	50.024 52.024
	WP	14.025	DESS PRE-MILLED ABUTMENT TI	DESS AST OSS PRE-MILLED TI	AST_OSS_WP E 61.025	50.025 52.025
Friadent XIVE®	NP	14.038	DESS PRE-MILLED ABUTMENT TI	DESS DEN XIV PRE-MILLED TI	DEN_XIV_NP E 61.038	50.038 52.038
	RP	14.039	DESS PRE-MILLED ABUTMENT TI	DESS DEN XIV PRE-MILLED TI	DEN_XIV_RP E 61.039	50.039 52.039
	WP	14.040	DESS PRE-MILLED ABUTMENT TI	DESS DEN XIV PRE-MILLED TI	DEN_XIV_WP E 61.040	50.040 52.040
Ankylos® C/X	All platforms, height GH 2,5	14.048 14.049	DESS PRE-MILLED ABUTMENT TI	DESS DEN ANK PRE-MILLED TI	DEN ANK E GH2 61.250	50.049 52.049
Osstem® TS	Mini Yellow	14.072	DESS PRE-MILLED ABUTMENT TI	DESS OSS TS PRE-MILLED TI	OSS_TS MINI_E 61.072	50.072 52.072
	STD Regular	14.073	DESS PRE-MILLED ABUTMENT TI	DESS OSS TS PRE-MILLED TI	OSS_TS STD_E 61.072	50.073 52.073



non-engaging



non-engaging



Scan Abutment
for desktop scanning
50.0XX



Scan Abutment
for intra-oral scanning
52.0XX

	Platform	Analogue	Category	System	Kit	Scan Abutments	
Camlog®	3.3	14.064	DESS PRE-MILLED ABUTMENT TI	DESS CAMLOG PRE-MILLED TI	CAM_3.3 E 61.064	52.064	
	3.8	14.065	DESS PRE-MILLED ABUTMENT TI	DESS CAMLOG PRE-MILLED TI	CAM_3.8 E 61.065	50.065	52.065
	4.3	14.066	DESS PRE-MILLED ABUTMENT TI	DESS CAMLOG PRE-MILLED TI	CAM_4.3 E 61.066	50.066	52.066
	5.0	14.067	DESS PRE-MILLED ABUTMENT TI	DESS CAMLOG PRE-MILLED TI	CAM_5.0 E 61.067	52.067	
Conelog®	NP	14.074	DESS PRE-MILLED ABUTMENT TI	DESS CONELOG PRE-MILLED TI	CONELOG NP E 61.074	52.074	
	RP	14.075	DESS PRE-MILLED ABUTMENT TI	DESS CONELOG PRE-MILLED TI	CONELOG RP E 61.075	52.075	
	WP	14.076	DESS PRE-MILLED ABUTMENT TI	DESS CONELOG PRE-MILLED TI	CONELOG WP E 61.076	52.076	
Biohorizns® External	RP/4.2	14.021	DESS PRE-MILLED ABUTMENT TI	DESS BH EXT PRE-MILLED TI	BH_EXTERNAL_RP E 61.021	50.021	52.021
	WP/5.2	14.022	DESS PRE-MILLED ABUTMENT TI	DESS BH EXT PRE-MILLED TI	BH_EXTERNAL_WP E 61.022	50.022 52.022	
MIS® Seven	SP	14.095	DESS PRE-MILLED ABUTMENT TI	DESS MIS SEV TI-BASE	MIS_SEV_SP_E 61.017	50.017	52.017
	WP	14.096	DESS PRE-MILLED ABUTMENT TI	DESS MIS SEV PRE-MILLED TI	MIS_SEV_WP_E 61.018	50.018	52.018
Megagen AntRidge®		14.057	DESS PRE-MILLED ABUTMENT TI	DESS MGN ANY PRE-MILLED TI	MGN ANY E 61.057	50.057	
Dentium Superline™ & Implantium	RP	14.090	DESS PRE-MILLED ABUTMENT TI	DESS DENTIUM PRE-MILLED TI	DENTIUM E 61.090	52.090	



non-engaging



non-engaging



Scan Abutment
for desktop scanning
50.0XX



Scan Abutment
for intra-oral scanning
52.0XX

**Table 1.5 DESS
Pre-milled
Abutment CoCr**

	Platform	Analogue	Category	System	Kit	Scan Abutments	
Nobel Active®	3.0	14.068	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB ACT PRE-MILLED CoCr	NOB_ACT_3.0 E 71.068	52.068	
	NP	14.041	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB ACT PRE-MILLED CoCr	NOB_ACT_NP E 71.041	50.041	52.041
	RP	14.042	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB ACT PRE-MILLED CoCr	NOB_ACT_RP E 71.042	50.042	52.042
	WP	14.069	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB ACT PRE-MILLED CoCr	NOB_ACT_WP E 71.069	52.069	
Nobel Replace Select™	NP	14.004	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB RPL PRE-MILLED CoCr	NOB_RPL_NP E 71.004	50.004	52.004
	RP	14.005	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB RPL PRE-MILLED CoCr	NOB_RPL_RP E 71.005	50.005	52.005
	WP	14.006	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB RPL PRE-MILLED CoCr	NOB_RPL_WP E 71.006	50.006	52.006
	6.0	14.097	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB RPL PRE-MILLED CoCr	NOB_RPL_NP 6.0 E 71.097	52.097	
Nobel Brånemark®	NP	14.001	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB BRK PRE-MILLED CoCr	NOB_BRK_NP E 71.001	50.001	52.001
	RP	14.002	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB BRK PRE-MILLED CoCr	NOB_BRK_RP E 71.002	50.002	52.002
	WP	14.003	DESS PRE-MILLED ABUTMENT CoCr	DESS NOB BRK PRE-MILLED CoCr	NOB_BRK_WP E 71.003	50.003 52.003	
3i Osseotite®	NP	14.011	DESS PRE-MILLED ABUTMENT CoCr	DESS 3i OSS PRE-MILLED CoCr	3i_OSS_NP E 71.011	50.011	52.011
	RP	14.012	DESS PRE-MILLED ABUTMENT CoCr	DESS 3i OSS PRE-MILLED CoCr	3i_OSS_RP E 71.012	50.012	52.012
	WP	14.013	DESS PRE-MILLED ABUTMENT CoCr	DESS 3i OSS PRE-MILLED CoCr	3i_OSS_WP E 71.013	50.013	52.013
3i Certain®	NP	14.014	DESS PRE-MILLED ABUTMENT CoCr	DESS 3i CER PRE-MILLED CoCr	3i_CER_NP E 71.014	50.014	52.014
	RP	14.015	DESS PRE-MILLED ABUTMENT CoCr	DESS 3i CER PRE-MILLED CoCr	3i_CER_RP E 71.015	50.015	52.015
	WP	14.016	DESS PRE-MILLED ABUTMENT CoCr	DESS 3i CER PRE-MILLED CoCr	3i_CER_WP E 71.016	50.016	52.016
Straumann® Soft Tissue Level	RN	14.009	DESS PRE-MILLED ABUTMENT CoCr	DESS STR STL PRE-MILLED CoCr	STR_STLV RN E 71.009	50.009 52.009	
	WN	14.045	DESS PRE-MILLED ABUTMENT CoCr	DESS STR STL PRE-MILLED CoCr	STR_STLV WN E 71.045	50.045 52.045	



non-engaging



non-engaging




Scan Abutment
for desktop scanning
50.0XX



Scan Abutment
for intra-oral scanning
52.0XX

	Platform	Analogue	Category	System	Kit	Scan Abutments	
Straumann® Bone Level	NC	14.043	DESS PRE-MILLED ABUTMENT CoCr	DESS STR BLV PRE-MILLED CoCr	STR_BLV_NC E 71.043	50.043	52.043
	RC	14.044	DESS PRE-MILLED ABUTMENT CoCr	DESS STR BLV PRE-MILLED CoCr	STR_BLV_RC E 71.044	50.044	52.044
Zimmer Screw-Vent®	NP	14.017	DESS PRE-MILLED ABUTMENT CoCr	DESS ZIM SCV PRE-MILLED CoCr	ZIM_SCV_NP E 71.017	50.017	52.017
	RP	14.018	DESS PRE-MILLED ABUTMENT CoCr	DESS ZIM SCV PRE-MILLED CoCr	ZIM_SCV_RP E 71.018	50.018	52.018
	WP	14.019	DESS PRE-MILLED ABUTMENT CoCr	DESS ZIM SCV PRE-MILLED CoCr	ZIM_SCV_WP E 71.019	50.019	52.019
Astra tech implant system™ EV	EV/3.6	14.059	DESS PRE-MILLED ABUTMENT CoCr	DESS AST EVO PRE-MILLED CoCr	AST_EVO 3.6_E 71.059	50.059	52.059
	EV/4.2	14.060	DESS PRE-MILLED ABUTMENT CoCr	DESS AST EVO PRE-MILLED CoCr	AST_EVO 4.2_E 71.060	50.060	52.060
	EV/4.8	14.061	DESS PRE-MILLED ABUTMENT CoCr	DESS AST EVO PRE-MILLED CoCr	AST_EVO 4.8_E 71.061	50.061	52.061
Astra tech implant Osseospeed™	RP	14.024	DESS PRE-MILLED ABUTMENT CoCr	DESS AST OSS PRE-MILLED CoCr	AST_OSS_RP E 71.024	50.024	52.024
	WP	14.025	DESS PRE-MILLED ABUTMENT CoCr	DESS AST OSS PRE-MILLED CoCr	AST_OSS_WP E 71.025	50.025	52.025
Friadent XIVE®	NP	14.038	DESS PRE-MILLED ABUTMENT CoCr	DESS DEN XIV PRE-MILLED CoCr	DEN_XIV_NP E 71.038	50.038	52.038
	RP	14.039	DESS PRE-MILLED ABUTMENT CoCr	DESS DEN XIV PRE-MILLED CoCr	DEN_XIV_RP E 71.039	50.039	52.039
	WP	14.040	DESS PRE-MILLED ABUTMENT CoCr	DESS DEN XIV PRE-MILLED CoCr	DEN_XIV_WP E 71.040	50.040	52.040
Ankylos® C/X	All platforms, height GH 2,5	14.048 14.049	DESS PRE-MILLED ABUTMENT CoCr	DESS DEN ANK PRE-MILLED CoCr	DEN ANK E GH2 71.250	50.049	52.049
Osstem® TS	Mini Yellow	14.072	DESS PRE-MILLED ABUTMENT CoCr	DESS OSS TS PRE-MILLED CoCr	OSS_TS MINI_E 71.072	50.072	52.072
	STD Regular	14.073	DESS PRE-MILLED ABUTMENT CoCr	DESS OSS TS PRE-MILLED CoCr	OSS_TS STD_E 71.072	50.073	52.073

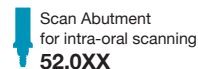

	Platform	Analogue	Category	System	Kit	Scan Abutments	
Camlog®	3.3	14.064	DESS PRE-MILLED ABUTMENT CoCr	DESS CAMLOG PRE-MILLED CoCr	CAM_3.3 E 71.064	52.064	
	3.8	14.065	DESS PRE-MILLED ABUTMENT CoCr	DESS CAMLOG PRE-MILLED CoCr	CAM_3.8 E 71.065	50.065	52.065
	4.3	14.066	DESS PRE-MILLED ABUTMENT CoCr	DESS CAMLOG PRE-MILLED CoCr	CAM_4.3 E 71.066	50.066	52.066
	5.0	14.067	DESS PRE-MILLED ABUTMENT CoCr	DESS CAMLOG PRE-MILLED CoCr	CAM_5.0 E 71.067	52.067	
Conelog®	NP	14.074	DESS PRE-MILLED ABUTMENT CoCr	DESS CONELOG PRE-MILLED CoCr	CONELOG NP E 71.074	52.074	
	RP	14.075	DESS PRE-MILLED ABUTMENT CoCr	DESS CONELOG PRE-MILLED CoCr	CONELOG RP E 71.075	52.075	
	WP	14.076	DESS PRE-MILLED ABUTMENT CoCr	DESS CONELOG PRE-MILLED CoCr	CONELOG WP E 71.076	52.076	
Biohorizns® External	RP/4.2	14.021	DESS PRE-MILLED ABUTMENT CoCr	DESS BH EXT PRE-MILLED CoCr	BH_EXTERNAL_RP E 671.021	50.021	52.021
	WP/5.2	14.022	DESS PRE-MILLED ABUTMENT CoCr	DESS BH EXT PRE-MILLED CoCr	BH_EXTERNAL_WP E 71.022	50.022 52.022	
MIS® Seven	SP	14.095	DESS PRE-MILLED ABUTMENT CoCr	DESS MIS SEV PRE-MILLED CoCr	MIS_SEV_SP_E 71.017	50.017	52.017
	WP	14.096	DESS PRE-MILLED ABUTMENT CoCr	DESS MIS SEV PRE-MILLED CoCr	MIS_SEV_WP_E 71.018	50.018	52.018
Megagen AntRidge®		14.057	DESS PRE-MILLED ABUTMENT CoCr	DESS MGN ANY PRE-MILLED CoCr	MGN ANY E 71.057	50.057	
Dentium Superline™ & Implantium	RP	14.090	DESS PRE-MILLED ABUTMENT CoCr	DESS DENTIUM PRE-MILLED CoCr	DENTIUM Ti base E 71.090	52.090 53.100	



non-engaging



non-engaging



Scan Abutment
for desktop scanning
50.0XX



Scan Abutment
for intra-oral scanning
52.0XX

All DESS products are manufactured according to ISO 9001:2015 and ISO 13485:2016 with CE marking for all products Class IIb

All DESS® AURUMBas® screws and screwdrivers are based on the standard ISO 10664

3i® and 3i Certain®, Miniplant®, MicroMiniplant™ and Osseotite® are registered trademarks of 3i Implant Innovations, USA

ASTRA TECH™, OsseoSpeed™, ANKYLOS® C/X, XiVE®, FRIADENT®, ASTRATECH implant system™ EV and UniAbutment

are registered trademarks of DENTSPLY SIRONA Implants

BioHorizons® is a registered trademark of BioHorizons, INC.

Brånemark System®, NobelReplace® and Replace Select™, NobelActive™, NobelReplace® CC, NobelSpeedy™ and Multi-Unit are

registered trademarks of Nobel Biocare AB, Sweden

ITF®, Straumann® and synOcta® are registered trademarks of Straumann Holding AG, Switzerland

MIS® is a registered trademark of MIS Implant Technologies Ltd.

Zimmer Screw-Vent® and SwissPlus® are registered trademarks of Zimmer Dental INC, USA

Dyna® is a registered trademark of Dyna Dental Engineering BV

Locator® is a registered trademark of ZEST IP Holdings

Camlog® and Conelog® are registered trademarks of CAMLOG Biotechnologies AG

Megagen® is a registered trademark of MEGAGEN IMPLANT CO. LTD

TORX® is a trademark of ACUMENT INTELLECTUAL PROPERTIES, LLC

Bego Semado® is a registered trademark of BEGO GmbH & Co. KG

Anthogyr® is a registered trademark of Anthogyr SAS

Medentis ICX is a registered trademark of medentis medical GmbH

Osstem® is a registered trademark of OSSTEM IMPLANT CO. LTD

Dentium® is a trademark of Dentium Co. Ltd.

GOBALD® is a trademark of Global D SASU.

Biotech Kontakt® is a trademark Biotech Dental SAS.

BTI® is a trademark BTI Biotechnology Institute.

DESS products can be used only by licensed dentists, physicians and labs.

Please refer ever to the instructions of use available on the website

www.dess-usa.com

or supplied with the product.

Some products may not be available in your country. Some of the products are cleared for US and Canada.

Please contact your local distributor for more details.

